Benchmark



Tube racks included with myBath 4L!



Technical Data and Ordering Information:

The myBath series is designed to guarantee superior accuracy and performance in a compact footprint. With four sizes to choose from, there is a myBath to fit every lab.

The 4, 8 and 12L models include easy to use controls allow the user to digitaly select and monitor settings. Adjustable in 0.1°C increments, the temperature is shown on a large, 3-digit display. The proprietary Quik-CAL[™] feature allows for simple, in-lab calibration.* High efficiency insulation maintains temperature stability and uniformity within the bath, while keeping the housing cool to the touch.

All myBath models feature seamless, stainless steel chambers for easy cleaning and long life The chambers resist rust and have no corners to collect dirt or contamination. The larger 8L and 12L units have drains for hassle-free water changes.


The 4L, 8L, and 12L models include a see-through cover that can be removed, or flipped back on a rear hinge to conserve bench space. The lid's gabled design will divert condensation to the front and back of the bath, away from samples.

Racks are available to accommodate 1.5/2.0ml, 15ml and 50 ml tubes within the 4L, 8L and 12L baths. The myBath 4L is supplied with all three racks.

*2L models use traditional calibration and include a flat, removable lid.

| | 2 Liter | 4 Liter | 8 Liter | 12 Liter | |
|-------------------|---------------------------------|------------------------|------------------------|------------------------|--|
| Catalog Number | B2000-2* | B2000-4* | B2000-8* | B2000-12* | |
| Temperature Range | Ambient + 5 to 100°C | | | | |
| Temp. Accuracy | ± 0.2°C at 37°C | | ±0.5°C at 37°C | | |
| Temp. Increments | 1.0°C | 0. | 0.1℃ | | |
| Temp. Uniformity | ± 0.2°C at 37°C | | | | |
| Exterior (WxDxH) | 6.5 x 8 x 6 in. | 7.7 x 13.9 x 9 in. | 11.2 x 17.5 x 12.4 in. | 11.2 x 17.5 x 16.2 in. | |
| Exterior (cm) | 16.5 x 20.3 x 15.2 cm. | 19.5 x 35.3 x 22.9 cm. | 28.4 x 44.5 x 31.5 cm. | 28.4 x 44.5 x 41.1 cm. | |
| Chamber (WxDxH) | 5.5 x 6 x 4 in. | 5.3 x 9.4 x 5.5 in. | 9.5 x 12.8 x 5.9 in. | 9.5 x 12.8 x 7.9 in. | |
| Chamber (cm) | 13.9 x 15.2 x 10.2 cm. | 13.5 x 23.9 x 13.9 cm. | 24.1 x 32.5 x 14.9 cm. | 24.1 x 32.5 x 20.1 cm. | |
| Weight | 4.5 lbs / 2.4 kg | 15 lbs / 6.8 kg | 17 lbs / 7.7 kg | 17 lbs / 7.7 kg | |
| Electrical | 115 or 230V, 50-60Hz, 250W-700W | | | | |

15V with US plug. For 230V (EU plug), please add (-E)

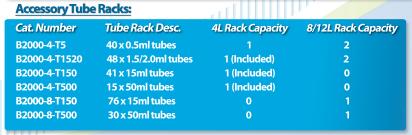



Benchmark

PO Box 709 Edison, NJ 08818 Phone: 908-769-5555 Fax: 732-313-7007 Web: www.BenchmarkScientific.com Email: Info@BenchmarkScientific.com